US MEDICA

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CONTENTS

Features.....	15
Important safety instructions.....	16
Operating instruction	18
Troubleshooting.....	22
Care & maintenance	23
Specifications.....	23

FEATURES

- Rolling foot massage
- For car, home and office
- Adjustable speed
- 2 rotation directions
- Timer: 15 minutes
- Includes AC and DC adapters for use at home or in a car
- Extra easy-to-clean cover and a storage bag included

IMPORTANT SAFETY INSTRUCTIONS

Please read this user's manual attentively.

- Do not put hands or fingers near the massage mechanism or rollers while this product is running.
- Touching the mechanism can result an injury caused by rollers moving.
- If the supply cord is damaged, it must be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid a hazard.
- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- If you feel any abnormal discomfort during the operation of this product, turn off the power immediately and do not use the product. Consult your physician.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into a water. Return this product to a service center for examination and repair.
- Always keep the power cord away from heated surfaces.
- Never drop or insert any object into any opening.
- When using this product after it has been in storage, check to see that the product is operating properly prior to use.
- To disconnect, turn all controls to the OFF position, then remove plug from outlet.
- Do not massage feet or calves that are swollen, inflamed or covered with skin eruptions.
- Connect this appliance to a properly grounded outlet only.
- Do not use this product if noise above normal sound is heard.
- Do not stand or put heavy objects on the massager. Use only while seated.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- It is not recommended to use the massager while driving.
- Plug the massager to 12V connector only after you start the car.

SAVE THESE INSTRUCTIONS

ATTENTION!

- **Do not** operate where aerosol (spray) products are being used or where oxygen is being administered.
- **Do not** operate under blanket or pillow.
- **Do not** stand on the massager.
- **Do not** pull by the AC cord.
- **Do not** carry this appliance by the power cord or use the cord as a handle.
- **Do not** reach for a product that has fallen into water. Unplug immediately.
- **Do not** use in the place with high humidity.


Persons who are undergoing medical treatment as well as the following persons should consult a physician before use:

- Persons who have injured their legs or feet in the past, persons who have numbness in their feet.
- Persons having a bone condition, such as osteoporosis.
- Persons having a circulatory disorder (heart disease, angiopathy, hypertension, etc.).
- Persons having a respiratory organ disorder.
- Persons using a body implanted medical electronic device.
- Persons having malignant tumors.
- Persons having symptoms of a circulatory disorder (thrombosis, serious aneurysm, chronic varicosis, etc.), dermatitis, or infectious skin diseases.
- Persons having a sensory disorder due to an advanced peripheral circulatory disorder, such as that resulting from diabetes.
- Persons having skin injuries.
- Persons having a fever (38°C or higher) due to disease.
- Persons who are pregnant or thought to be pregnant or who are menstruating.
- Persons who feel physical discomfort and require rest.
- Persons wishing to use the unit for rehabilitation purposes.
- Persons other than the above who do not feel physically normal.

Consult your physician
before use!

OPERATING INSTRUCTION

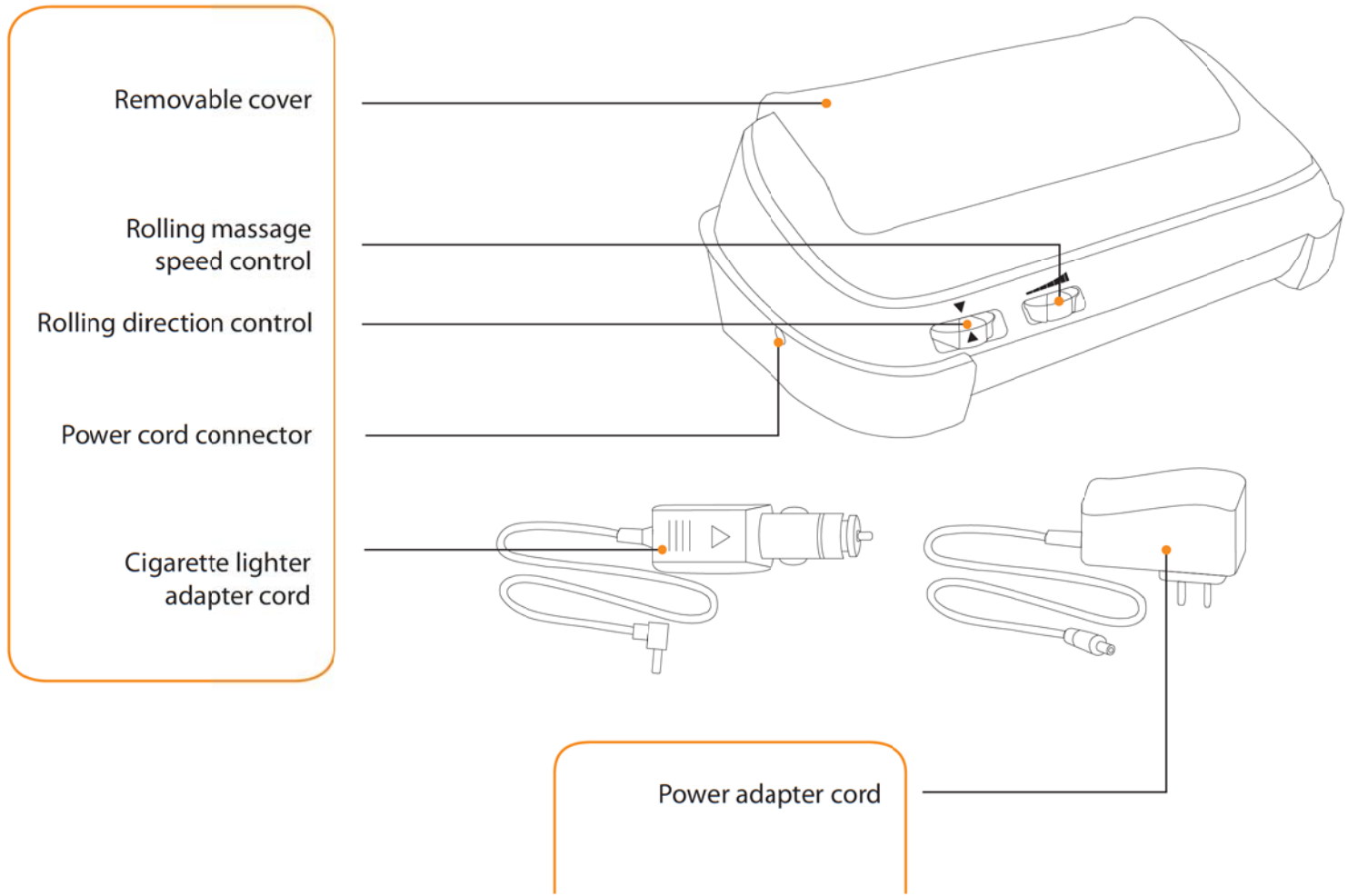IN A CAR

AT HOME

Place the massager on a flat surface.



Plug one end of the cigarette lighter adapter cord into the power cord connector on the massager. Start the car and then plug the other end into the 12V cigarette lighter connector in the car.

Plug one end of the power adapter cord into the power cord connector on the massager. Plug the other end into the 220V socket.


Adjust the direction control to choose one of the rotation directions: clockwise «▲» or counter clockwise «▼». The massage rollers will start movement. Adjust the rollers' rotation speed using the speed control if necessary. In case you need to stop the massage immediately, set the direction control to the «■» position.

Cover removing

1. Before taking off the removable cover, make sure that the massager is unplugged. Using a sharp object pry up from both sides the metal pivot, located on the back side of the massager and remove it off the holding bolts (see the picture to the right).

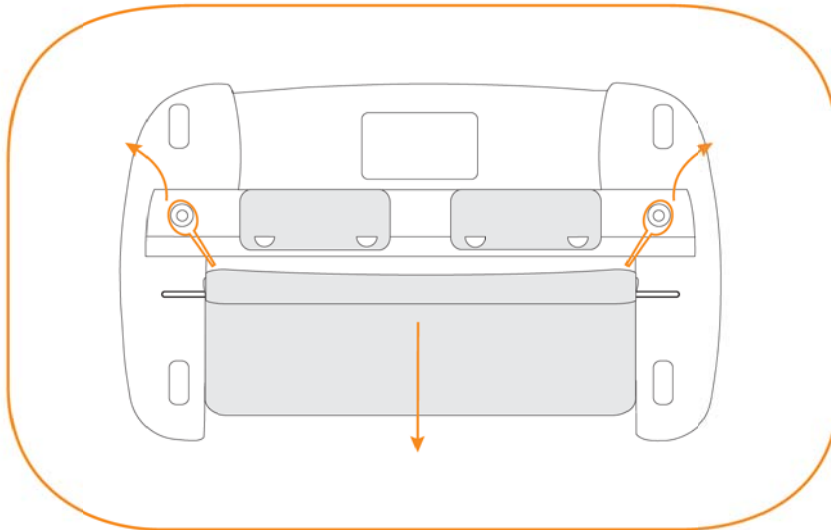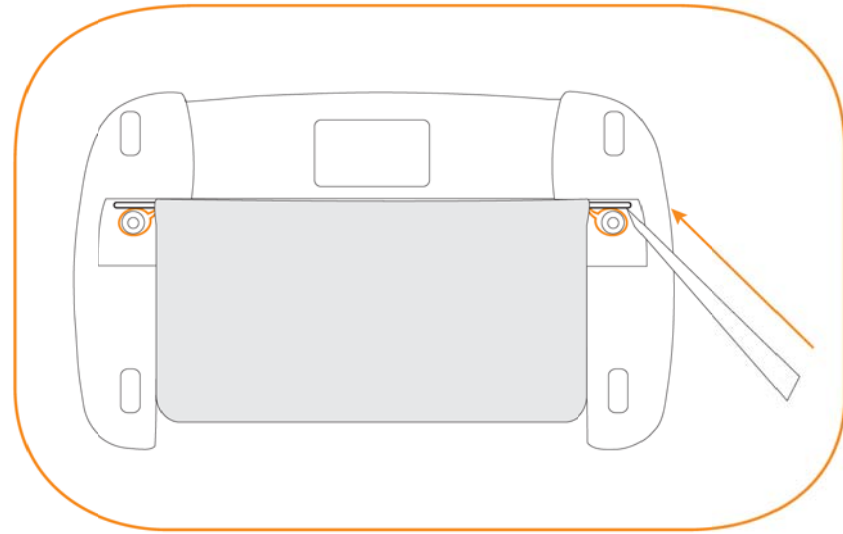
2. Take the loops off the holding bolts and throw the cover, as shown on the picture to the left.

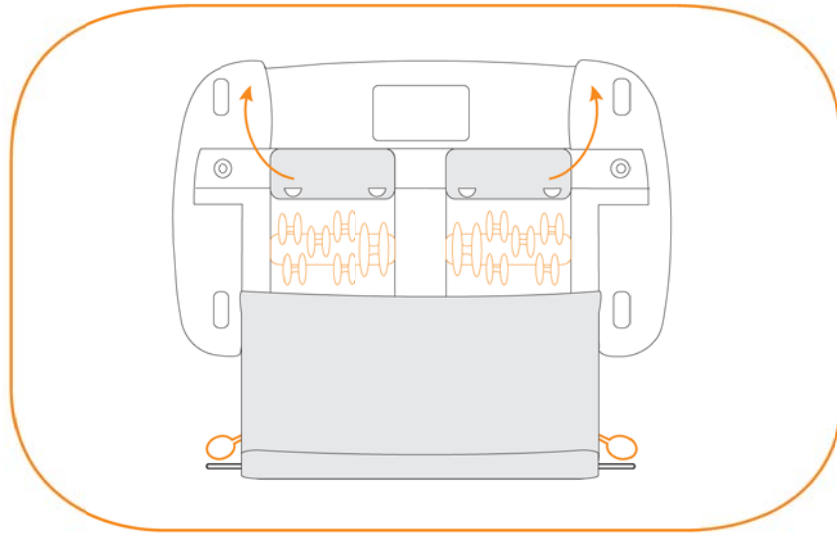3. Take 4 remaining loops off the plastic shafts, located in front of the massage mechanism (see the picture to the left).

4. Pull the cover to remove it, as shown on the picture to the right.

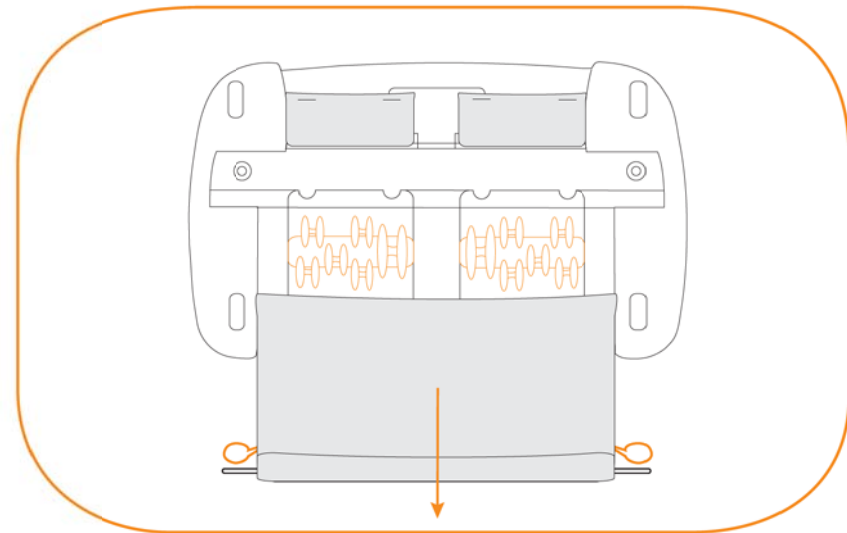TROUBLESHOOTING

PROBLEM	CAUSE	SOLUTION
Massage stopped	The preset maximum period of operation has elapsed	Let the device cool down for at least 15 min before using again. Switch the device on again
	The connection to the mains has been interrupted	Check the plug connection and the contact of the main plug with the main supply
	Rolling direction control is set to the «■» position	Set the rolling direction control to the «▼» or «▲» position
Loud noise occurs during the massage	The noise can be caused by the massage mechanisms	The noise is not a trouble of the massager

CARE & MAINTENANCE

Storage

- Do not wrap the power cord around the device before storing.
- Store the foot massager in its storage bag and in a safe, dry and cool place.
- Store it out of the reach of children.
- Avoid contact with sharp edges or pointed objects as this might damage the fabric of the cover.

Cleaning

- Unplug the device.
- Use a soft, damp cloth to clean the device.
- Use a soft brush to clean dirty connections and contacts.
- Dry the device with a soft cloth to minimise the potential water damage.
- Wait until the device is completely dry before using it again.

SPECIFICATIONS

Name of the product	Foot massager
Model	OMEGA
Voltage	Entrance: 220V, 12V
Power	19 W
Dimensions (Height/Length/Width)	4,3" x 10" x 5,6"
Timer	15 minutes

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US Medica LLC., 1201 N.Orange St, Wilmington, DE 19801, USA